

ZAP70 antibody [MSVA-870R] HistoMAX™

Cat. No. GTX04441

Host	Rabbit
Clonality	Monoclonal
Isotype	IgG
Applications	IHC-P
Reactivity	Human

Package

500 µl, 100 µl

PRODUCT

This antibody was validated on 76 different Normal Tissues by IHC-P.

Summary

[Go to Normal Tissue Gallery](#)[Go to Cancer Tissue Gallery](#)

Applications

Application Note

*Optimal dilutions/concentrations should be determined by the researcher.

Suggested dilution	Recommended dilution
IHC-P	1:75-1:150

Note : Manual staining : Heat-induced antigen retrieval for 5 minutes in an autoclave at 121°C in pH 7.8 Tris-EDTA-based Target Retrieval Solution buffer.

Not tested in other applications.

Product Note **Highly recommended for IHC-P in human tissues.**

Properties

Form	Liquid
Buffer	PBS, 0.05% BSA (Please contact us for PBS only format)
Preservative	0.05% sodium azide
Storage	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
Concentration	0.2 mg/ml (Please refer to the vial label for the specific concentration.)
Immunogen	Recombinant fragment (around aa 247-382) of human ZAP70 protein
Purification	Protein A/G purified
Conjugation	Unconjugated



For full product information, images and publications, please visit our [website](#).


Date 2026 / 10 / 10 Page 1 of 2

For laboratory research use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.

Note


Purchasers shall not, and agree not to enable third parties to, analyze, copy, reverse engineer or otherwise attempt to determine the structure or sequence of the product.

DATA IMAGES

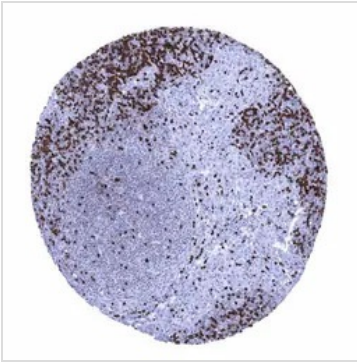

GTx04441 IHC-P Image

IHC-P analysis of human lymph node tissue using GTx04441 ZAP70 antibody [MSVA-870R] HistoMAX™. In lymph nodes ZAP70 positive lymphocytes predominate in the interfollicular zone but are much less numerous in germinal centres.

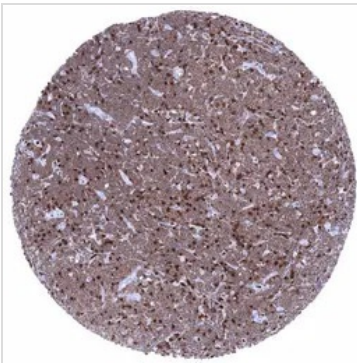

GTx04441 IHC-P Image

IHC-P analysis of human B-cell chronic lymphocytic leukemia (B-CLL) tissue using GTx04441 ZAP70 antibody [MSVA-870R] HistoMAX™.

B-CLL with strong ZAP70 staining of tumor cells and even stronger staining in interspersed normal T lymphocytes.



GTx04441 IHC-P Image

IHC-P analysis of human B-cell chronic lymphocytic leukemia (B-CLL) tissue using GTx04441 ZAP70 antibody [MSVA-870R] HistoMAX™.

ZAP70 negative B-CLL containing a significant number of ZAP70 positive T lymphocytes.



For full product information, images and publications, please visit our [website](#).