

# Nucleolin antibody [MSVA-623R] HistoMAX™

**Cat. No. GTX04456**

<b>Host</b>	Rabbit
<b>Clonality</b>	Monoclonal
<b>Isotype</b>	IgG
<b>Applications</b>	IHC-P
<b>Reactivity</b>	Human

## Package

500 µl, 100 µl

## PRODUCT

This antibody was validated on 76 different Normal Tissues by IHC-P.

## Summary

[Go to Normal Tissue Gallery](#)

[Go to Cancer Tissue Gallery](#)

## Applications

### Application Note

\*Optimal dilutions/concentrations should be determined by the researcher.

### Suggested dilution Recommended dilution

IHC-P 1:100-1:200

**Note : Manual staining : Heat-induced antigen retrieval for 5 minutes in an autoclave at 121°C in pH 7.8 Tris-EDTA-based Target Retrieval Solution buffer.**

**This antibody has been validated by Autostainer including Agilent/Dako - Autostainer Link 48, Bond RX-Lecia, and Roche-Ventana Discovery ULTRA. The detailed steps by Autostainer please refer to the [protocol](#).**

Not tested in other applications.

## Product Note

**Highly recommended for IHC-P in human tissues.** Autostainer protocol information available

## Properties

<b>Form</b>	Liquid
<b>Buffer</b>	PBS, 0.05% BSA (Please contact us for PBS only format)
<b>Preservative</b>	0.05% sodium azide
<b>Storage</b>	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
<b>Concentration</b>	0.2 mg/ml (Please refer to the vial label for the specific concentration.)
<b>Immunogen</b>	Recombinant full-length human NCL protein
<b>Purification</b>	Protein A/G purified



For full product information, images and publications, please visit our [website](#).

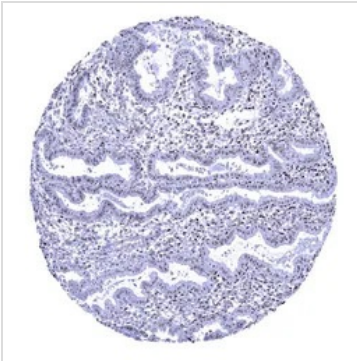
**Conjugation**

Unconjugated

**Note**

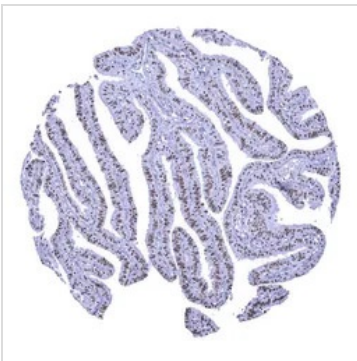
For laboratory research use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.

Purchasers shall not, and agree not to enable third parties to, analyze, copy, reverse engineer or otherwise attempt to determine the structure or sequence of the product.

**DATA IMAGES**

**GTx04456 IHC-P Image**

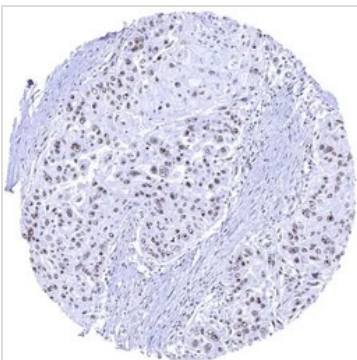
IHC-P analysis of human endometrium tissue using GTx04456 Nucleolin antibody [MSVA-623R] HistoMAX™.

In this sample of endometrium nucleolin staining is particularly prominent in the nucleoli of epithelial and stromal cells.


**GTx04456 IHC-P Image**

IHC-P analysis of human fallopian tube mucosa tissue using GTx04456 Nucleolin antibody [MSVA-623R] HistoMAX™.

The entire nuclei of epithelial cells are nucleolin positive but the immunostaining is most distinct in the nucleolus.


**GTx04456 IHC-P Image**

IHC-P analysis of human esophageal squamous cell carcinoma (ESCC) tissue using GTx04456 Nucleolin antibody [MSVA-623R] HistoMAX™.

Distinct nucleolin immunostaining in the nucleoli of an esophageal squamous cell carcinoma.



For full product information, images and publications, please visit our [website](https://www.genetex.com).