

Periostin antibody [MSVA-649M] HistoMAX™

Cat. No. GTX04460

Host	Mouse
Clonality	Monoclonal
Isotype	IgG2b
Applications	IHC-P
Reactivity	Human

Package

500 µl, 100 µl

PRODUCT

This antibody was validated on 76 different Normal Tissues by IHC-P.

Summary

[Go to Normal Tissue Gallery](#)

[Go to Cancer Tissue Gallery](#)

Applications

Application Note

*Optimal dilutions/concentrations should be determined by the researcher.

Suggested dilution	Recommended dilution
IHC-P	1:100-1:200

Note : Manual staining : Heat-induced antigen retrieval for 5 minutes in an autoclave at 121°C in pH 7.8 Tris-EDTA-based Target Retrieval Solution buffer.

Not tested in other applications.

Product Note

Highly recommended for IHC-P in human tissues.

Properties

Form	Liquid
Buffer	PBS, 0.05% BSA (Please contact us for PBS only format)
Preservative	0.05% sodium azide
Storage	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
Concentration	0.2 mg/ml (Please refer to the vial label for the specific concentration.)
Immunogen	Recombinant fragment (around aa193-326) of human POSTN protein
Purification	Protein A/G purified
Conjugation	Unconjugated



For full product information, images and publications, please visit our [website](#).


Date 2026 / 10 / 09 Page 1 of 2

For laboratory research use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.

Note


Purchasers shall not, and agree not to enable third parties to, analyze, copy, reverse engineer or otherwise attempt to determine the structure or sequence of the product.

DATA IMAGES

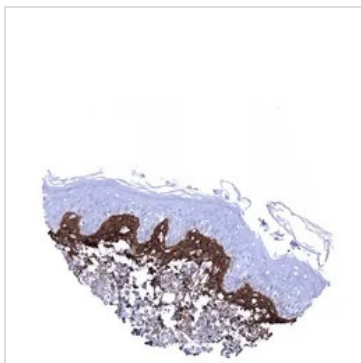


GTx04460 IHC-P Image

IHC-P analysis of human skin tissue using GTx04460 Periostin antibody [MSVA-649M] HistoMAX™. Periostin is most regularly seen in a subepithelial band of variable thickness in the skin.

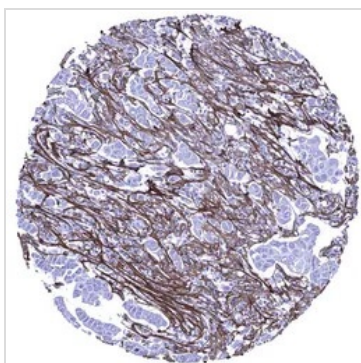

GTx04460 IHC-P Image

IHC-P analysis of human invasive lobular breast cancer (ILC) tissue using GTx04460 Periostin antibody [MSVA-649M] HistoMAX™. Abundant periostin positive stroma in an invasive lobular breast cancer.

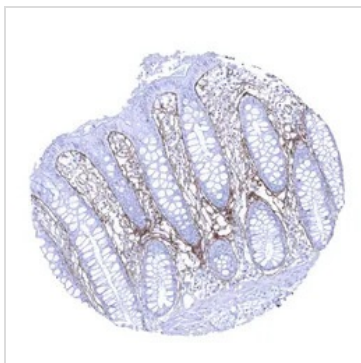

GTx04460 IHC-P Image

IHC-P analysis of human mucosa from rectum tissue using GTx04460 Periostin antibody [MSVA-649M] HistoMAX™. A periostin positive subepithelial band is regularly seen in the colorectal mucosa. Additional periostin deposits occur in the lamina propria.



For full product information, images and publications, please visit our [website](https://www.genetex.com).