

TLE1 antibody [MSVA-856R] HistoMAX™

Cat. No. GTX04966

Host	Rabbit
Clonality	Monoclonal
Isotype	IgG
Applications	IHC-P
Reactivity	Human

Package

500 µl, 100 µl

PRODUCT

Summary

[Go to Normal Tissue Gallery](#)

[Go to Cancer Tissue Gallery](#)

Applications

Application Note

*Optimal dilutions/concentrations should be determined by the researcher.

Suggested dilution

Recommended dilution

IHC-P

1:100 - 1:200

Note : Heat-induced antigen retrieval for 5 minutes in an autoclave at 121°C in pH 7.8 Tris-EDTA-based Target Retrieval Solution buffer.

Not tested in other applications.

Product Note

Highly recommended for IHC-P in human tissues.

Properties

Form	Liquid
Buffer	10mM PBS, BSA (Please contact us for sodium azide-free or sodium azide and BSA-free format)
Preservative	sodium azide
Storage	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
Concentration	66.66 µg/ml (Please refer to the vial label for the specific concentration.)
Immunogen	Recombinant human TLE1 fragment
Purification	Protein A/G purified
Conjugation	Unconjugated



For full product information, images and publications, please visit our [website](#).

Date 2026 / 10 / 07 Page 1 of 2

For laboratory research use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.

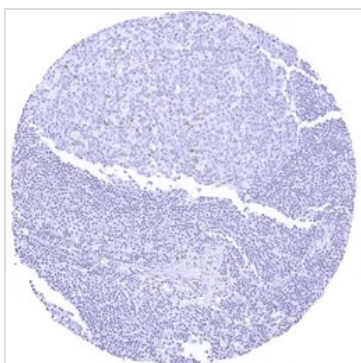
Note

Purchasers shall not, and agree not to enable third parties to, analyze, copy, reverse engineer or otherwise attempt to determine the structure or sequence of the product.

DATA IMAGES

GTX04966 IHC-P Image

IHC-P analysis of human endometrium tissue using GTX04966 TLE1 antibody [MSVA-856R] HistoMAX™. Endometrium with moderate to strong TLE1 staining of epithelial and stromal cells.


GTX04966 IHC-P Image

IHC-P analysis of human lymph node tissue using GTX04966 TLE1 antibody [MSVA-856R] HistoMAX™. Lymph node with weak to moderate TLE1 staining in a subset of lymphocytes predominantly in germinal centres.


GTX04966 IHC-P Image

IHC-P analysis of human synovial sarcoma tissue using GTX04966 TLE1 antibody [MSVA-856R] HistoMAX™. Synovial sarcoma with strong TLE1 positivity of all tumor cells.



For full product information, images and publications, please visit our [website](#).