

Influenza B virus Nucleoprotein (B/Sydney/3/2004), DDDDK Tag

Cat. No. GTX135867-pro

Application ELISA, Lateral Flow, Sandwich ELISA**Species** Influenza B virus

Reference (2)

★★★★★ Review (2)

Package
100 µg

APPLICATION

Application Note

Recommended antibody pairs for sandwich ELISA:

Capture : GTX636100 / GTX636099 / GTX636194, Detection : GTX636099 / GTX636100 / GTX636194 ; Capture: GTX629882 / GTX629883, Detection : GTX128539, or Capture : GTX128539, Detection : GTX629882

Please notice that detection antibodies need to be conjugated to HRP when paired with capture antibodies. Please contact us for custom HRP-conjugated antibody.

Recommended antibody pairs for Lateral Flow:

Capture : GTX636099 / GTX636100 / GTX636194, Detection: GTX636100 / GTX636099

Please notice that the detection antibodies need to be conjugated to Gold to function when paired with capture antibodies. Please contact us for custom Gold-conjugated antibody.

PROPERTIES

Form Liquid**Buffer** 90mM Tris, 1M Arginine**Preservative** No preservative**Storage** Store as concentrated solution. Centrifuge briefly prior to opening vial. Aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles. For long-term storage after reconstitution, aliquot and store at -70°C or below. Do not vortex.**Concentration** 0.8 mg/ml (Please refer to the vial label for the specific concentration.)**Region/Sequence** DDDDK tagged full length Nucleoprotein protein of Influenza B virus (B/Sydney/3/2004) (AGX16066)**Expression System** HEK293 cells**Purity** >90%**Conjugation** Unconjugated

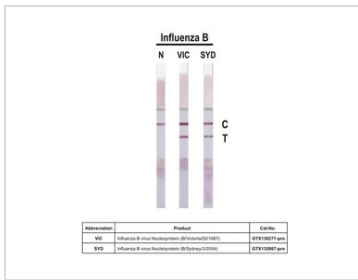
Note

For laboratory research use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.

Purchasers shall not, and agree not to enable third parties to, analyze, copy, reverse engineer or otherwise attempt to determine the structure or sequence of the product.



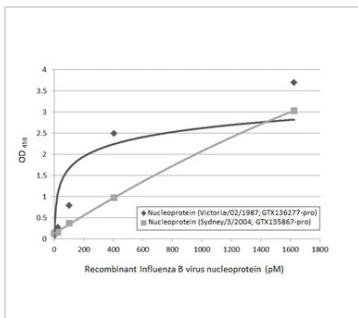
For full product information, images and publications, please visit our [website](#).

DATA IMAGES

GTX135867-pro Lateral Flow Image

Detection of recombinant influenza B virus nucleoproteins of the indicated strains by lateral flow assay using the recombinant rabbit monoclonal antibody pair.

Capture: Influenza B virus Nucleoprotein antibody [HL1069] (GTX636100)

Detection: Influenza B virus Nucleoprotein antibody [HL1068] (GTX636099)

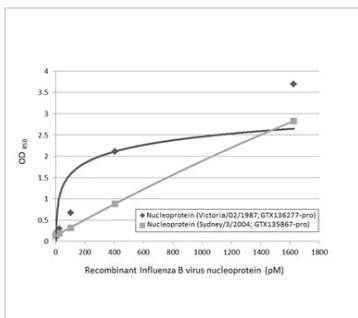

GTX135867-pro ELISA Image

Sandwich ELISA detection of recombinant nucleoproteins (NP) derived from different strains of Influenza B virus (i.e., B/Sydney/3/2004; B/Victoria/02/1987) using antibodies as below.

Capture: Influenza B virus Nucleoprotein antibody [HL1068] (GTX636099) (5 µg/mL)

Detection: HRP-conjugated Influenza B virus Nucleoprotein antibody [HL1069] (GTX636100) (1 µg/mL).

Please notice that GTX636100 needs to be conjugated to HRP to function as the detection antibody when paired with GTX636099. Please contact us for custom HRP-conjugated antibody.

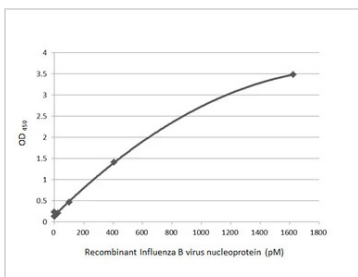

GTX135867-pro ELISA Image

Sandwich ELISA detection of recombinant nucleoproteins (NP) derived from different strains of Influenza B virus (i.e., B/Sydney/3/2004; B/Victoria/02/1987) using antibodies as below.

Capture: Influenza B virus Nucleoprotein antibody [HL1073] (GTX636194) (5 µg/mL)

Detection: HRP-conjugated Influenza B virus Nucleoprotein antibody [HL1068] (GTX636099) (1 µg/mL).

Please notice that GTX636099 needs to be conjugated to HRP to function as the detection antibody when paired with GTX636194. Please contact us for custom HRP-conjugated antibody.


GTX135867-pro ELISA Image

Sandwich ELISA detection of recombinant Influenza B virus Nucleoprotein (B/Sydney/3/2004), DDDDK Tag (GTX135867-pro) using antibodies as below.

Capture: Influenza B virus Nucleoprotein antibody [GT371] (GTX629882) (5 µg/mL)

Detection: Influenza B virus Nucleoprotein antibody (GTX128539) (1 µg/mL).



For full product information, images and publications, please visit our [website](https://www.genetex.com).