

NF-L antibody

Cat. No. GTX30190

| | |
|--------------|---------------------------------------|
| Host | Rabbit |
| Clonality | Polyclonal |
| Isotype | IgG |
| Applications | WB, ICC/IF, IHC-P |
| Reactivity | Human, Mouse, Rat, Bovine, Pig, Avian |

Package

50 µl

Applications

Application Note

*Optimal dilutions/concentrations should be determined by the researcher.

| Suggested dilution | Recommended dilution |
|--------------------|----------------------|
| WB | 1:20000 |
| ICC/IF | 1:500 - 1:5000 |
| IHC-P | 1:5000 |

Not tested in other applications.

Product Note Specifically recognizes the light neurofilament subunit (~68 kDa).

Properties

| | |
|--------------|--|
| Form | Liquid |
| Buffer | Serum |
| Preservative | No preservatives |
| Storage | Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles. |
| Immunogen | Native Neurofilament Light protein purified from pig spinal cord [UniProt# P02547] |
| Purification | Unpurified |
| Conjugation | Unconjugated |

Note


For laboratory research use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.

Purchasers shall not, and agree not to enable third parties to, analyze, copy, reverse engineer or otherwise attempt to determine the structure or sequence of the product.



For full product information, images and publications, please visit our [website](#).


DATA IMAGES

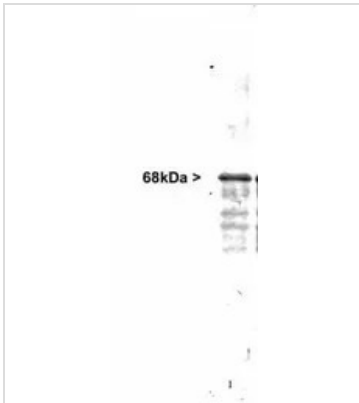

GTX30190 WB Image

WB analysis of rat cerebellum tissue using GTX30190 NF-L antibody.

Dilution : 1:20000

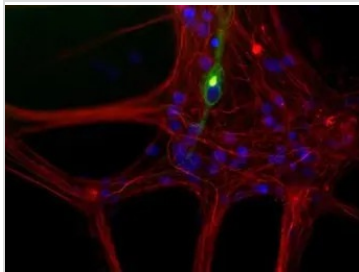

GTX30190 ICC/IF Image

ICC/IF analysis of mixed neuron/glia cell cultures from newborn rat brain using GTX30190 NF-L antibody (red), co-stained with anti-peripherin antibody (green) and DAPI . A class of large neurons, like the one in the middle of this image, contain peripherin, while the majority of neurons and their processes contain NF-L and not peripherin. Interestingly, the peripherin positive cells often contain a cytoplasmic inclusion next to the nucleus which stains for both peripherin and NF-L, and so appears golden in this kind of image.



For full product information, images and publications, please visit our [website](#).