

Survivin antibody [Mix]

Cat. No. GTX34059

| | |
|--------------|-------------------|
| Host | Mouse |
| Clonality | Multiclonal |
| Isotype | IgG1 |
| Applications | WB, IHC-P |
| Reactivity | Human, Mouse, Rat |

Package
100 µl

Applications

Application Note

*Optimal dilutions/concentrations should be determined by the researcher.

| Suggested dilution | Recommended dilution |
|--------------------|----------------------|
| WB | 1:1000-1:2000 |
| IHC-P | 1:200-1:500 |

Not tested in other applications.

Calculated MW 16 kDa. ([Note](#))

Properties

| | |
|---------------|--|
| Form | Liquid |
| Buffer | PBS, 0.5% BSA, 50% Glycerol |
| Preservative | 0.02% Sodium azide |
| Storage | Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles. |
| Concentration | Batch dependent (Please refer to the vial label for the specific concentration.) |
| Immunogen | Recombinant Protein |
| Purification | Purified by antigen-affinity chromatography |
| Conjugation | Unconjugated |

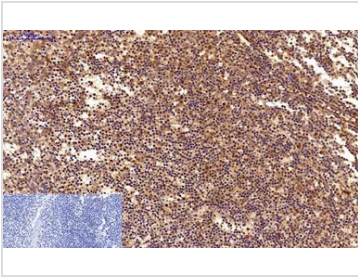
Note

For laboratory research use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.

Purchasers shall not, and agree not to enable third parties to, analyze, copy, reverse engineer or otherwise attempt to determine the structure or sequence of the product.



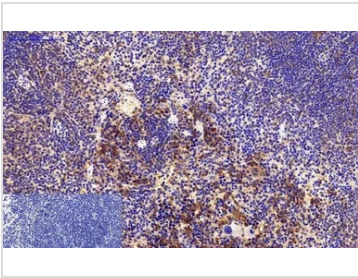
For full product information, images and publications, please visit our [website](#).

DATA IMAGES

GTX34059 IHC-P Image

IHC-P analysis of human tonsil tissue using GTX34059 Survivin antibody [Mix]. Negative control (the lower left coner) was secondary antibody only.

Antigen retrieval : Sodium citrate pH6.0 was used for antibody retrieval (>98°C, 20min)

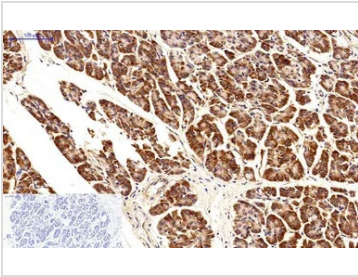
Dilution : 1:200


GTX34059 IHC-P Image

IHC-P analysis of mouse spleen tissue using GTX34059 Survivin antibody [Mix]. Negative control (the lower left coner) was secondary antibody only.

Antigen retrieval : Sodium citrate pH6.0 was used for antibody retrieval (>98°C, 20min)

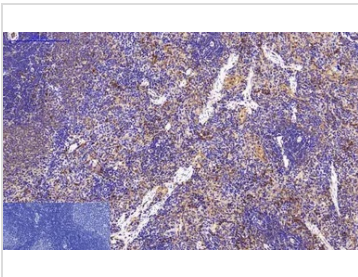
Dilution : 1:200


GTX34059 IHC-P Image

IHC-P analysis of human stomach cancer tissue using GTX34059 Survivin antibody [Mix]. Negative control (the lower left coner) was secondary antibody only.

Antigen retrieval : Sodium citrate pH6.0 was used for antibody retrieval (>98°C, 20min)

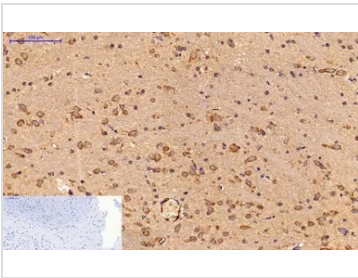
Dilution : 1:200


GTX34059 IHC-P Image

IHC-P analysis of rat spleen tissue using GTX34059 Survivin antibody [Mix]. Negative control (the lower left coner) was secondary antibody only.

Antigen retrieval : Sodium citrate pH6.0 was used for antibody retrieval (>98°C, 20min)

Dilution : 1:200


GTX34059 IHC-P Image

IHC-P analysis of rat brain tissue using GTX34059 Survivin antibody [Mix]. Negative control (the lower left coner) was secondary antibody only.

Antigen retrieval : Sodium citrate pH6.0 was used for antibody retrieval (>98°C, 20min)

Dilution : 1:200



For full product information, images and publications, please visit our [website](https://www.genetex.com).