

SARS-CoV-2 (COVID-19) Spike RBD ELISA Pair [HL1014 / HL1003]

Cat. No. GTX500046

Application	ELISA, Sandwich ELISA	Package 1 pair
Reactivity	SARS Coronavirus 2	·

PRODUCT

Product Content

Cat No.	Туре	Host	Conjugation	Package
GTX635792-01	Detection Antibody	Rabbit	HRP	50μΙ
GTX635807	Capture Antibody	Rabbit	Unconjugated	100μΙ

APPLICATION

Application Note

*Optimal dilutions/concentrations should be determined by the researcher.

Suggested dilution	Recommended dilution
ELISA	Assay dependent
Sandwich ELISA	Assay dependent

Not tested in other applications.

Product Note

This ELISA pair detects SARS-CoV-2 Spike protein, but does not cross-react with SARS-CoV or MERS-CoV spike proteins based on our internal testing. This ELISA pair is able to detect multiple SARS-CoV-2 VOCs, including Omicron variants (BA.1 and BA.2).

PROPERTIES	
Form	Liquid
Storage	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
Note	For <i>In vitro</i> laboratory use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.




For full product information, images and publications, please visit our <u>website</u>.

Date 2024 / 03 / 28 Page 1 of 2



DATA IMAGES

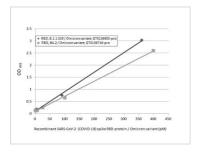


GTX500046 ELISA Image

Sandwich ELISA detection of recombinant Spike RBD Protein(s) derived from different strains of SARS-CoV-2 virus (ie., B.1.1.529 Omicron variant; BA.2 Omicron variant) using antibodies as below.

Capture: SARS-CoV-2 (COVID-19) Spike RBD antibody [HL1014] (GTX635807) (5 μg/mL)

Detection: SARS-CoV-2 (COVID-19) Spike RBD antibody [HL1003] (HRP) (GTX635792-01) (1 μg/mL)

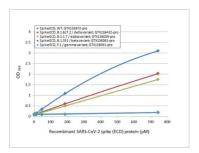


GTX500046 ELISA Image

Sandwich ELISA detection of recombinant Spike ECD protein(s) derived from different strains of SARS-CoV-2 virus (ie., Wild type; B.1.1.7 alpha variant; B.1.351 beta variant; P.1 gamma variant; B1.617.2 delta variant) using antibodies as below.


 $\textbf{Capture:} \ \mathsf{SARS-CoV-2} \ (\mathsf{COVID-19}) \ \mathsf{Spike} \ \mathsf{RBD} \ \mathsf{antibody} \ [\mathsf{HL1014}] \ (\mathsf{GTX635807}) \ (5 \ \mu \mathsf{g/mL})$

Detection: SARS-CoV-2 (COVID-19) Spike RBD antibody [HL1003] (HRP) (GTX635792-01) (1 μg/mL)

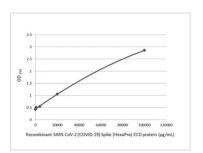


GTX500046 ELISA Image

Sandwich ELISA detection of recombinant SARS-CoV-2 (COVID-19) Spike (HexaPro) ECD, His tag (GTX135972-pro) using SARS-CoV-2 (COVID-19) Spike RBD antibody [HL1014] (GTX635807) as capture antibody at concentration of 5 μ g/mL and SARS-CoV-2 (COVID-19) Spike RBD antibody [HL1003] (HRP) (GTX635792-01) as detection antibody at concentration of 1 μ g/mL

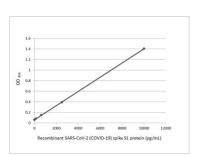

GTX500046 ELISA Image

Sandwich ELISA detection of recombinant SARS-CoV-2 (COVID-19) Spike S1 protein, His tag (active) protein (GTX135817-pro) using SARS-CoV-2 (COVID-19) Spike RBD antibody [HL1014] (GTX635807) as capture antibody at concentration of 5 μ g/mL and HRP-conjugated SARS-CoV-2 (COVID-19) Spike RBD antibody [HL1003] (GTX635792-01) as detection antibody at concentration of 1 μ g/mL

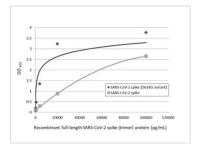

GTX500046 ELISA Image

Sandwich ELISA detection of recombinant SARS-CoV-2 (COVID-19) Spike (D614G variant) protein, His tag (active) (GTX02575-pro), and SARS-CoV-2 spike (trimer) protein using SARS-CoV-2 (COVID-19) Spike RBD antibody [HL1014] (GTX635807) as capture antibody at concentration of 5 μ g/mL and SARS-CoV-2 (COVID-19) Spike RBD antibody [HL1003] (HRP) (GTX635792-01) as detection antibody at concentration of 1 μ g/mL.



For full product information, images and publications, please visit our <u>website</u>.

Date 2024 / 03 / 28 Page 2 of 2