

KDM6A antibody [HMV311] HistoMAX™

Cat. No. GTX639935

Host	Rabbit
Clonality	Monoclonal
Isotype	IgG
Applications	IHC-P
Reactivity	Human

References (2)

Package

500 µl, 100 µl

PRODUCT

This antibody was validated on 76 different Normal Tissues by IHC-P.

Summary

[Go to Normal Tissue Gallery](#)[Go to Cancer Tissue Gallery](#)

Applications

Application Note

*Optimal dilutions/concentrations should be determined by the researcher.

Suggested dilution	Recommended dilution
--------------------	----------------------

IHC-P	1:100-1:200
-------	-------------

Note : Heat-induced antigen retrieval for 5 minutes in an autoclave at 121°C in Antigen retrieval buffer, pH 7.8 (GTX443621).

Not tested in other applications.

Product Note **Highly recommended for IHC-P in human tissues.**

Properties

Form	Liquid
Buffer	PBS, 2% BSA
Preservative	0.025% ProClin 300
Storage	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
Concentration	0.77 mg/ml
Immunogen	Recombinant protein encompassing a sequence within the center region of human KDM6A. The exact sequence is proprietary.
Purification	Affinity purified by Protein A.
Conjugation	Unconjugated



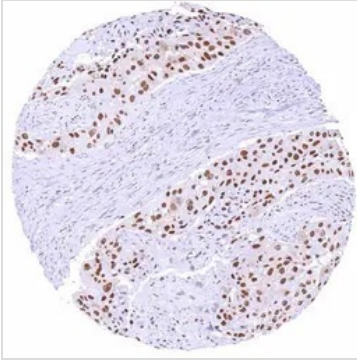
For full product information, images and publications, please visit our [website](#).

Date 2026 / 10 / 10 Page 1 of 2

For laboratory research use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.

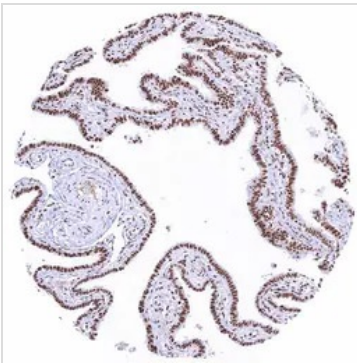
Note

Purchasers shall not, and agree not to enable third parties to, analyze, copy, reverse engineer or otherwise attempt to determine the structure or sequence of the product.

DATA IMAGES

GTx639935 IHC-P Image

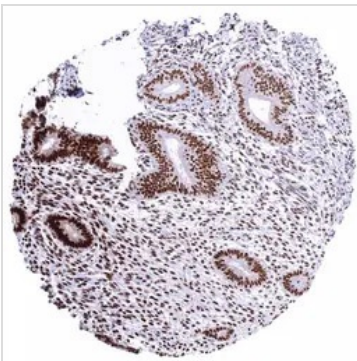
IHC-P analysis of human squamous cell carcinoma (SCC) from penile tissue using GTx639935 KDM6A antibody [HmV311] HistoMAX™.

Squamous cell carcinoma with strong KDM6A staining of tumor cells.


GTx639935 IHC-P Image

IHC-P analysis of human fallopian tube mucosa tissue using GTx639935 KDM6A antibody [HmV311] HistoMAX™.

Strong nuclear KDM6A staining of all epithelial cells.


GTx639935 IHC-P Image

IHC-P analysis of human endometrium (proliferation) tissue using GTx639935 KDM6A antibody [HmV311] HistoMAX™.

Strong nuclear KDM6A staining of epithelial and stromal cells.



For full product information, images and publications, please visit our [website](https://www.genetex.com).