

NFIX antibody [HMV329] HistoMAX™

Cat. No. GTX640561

Host	Rabbit
Clonality	Monoclonal
Isotype	IgG
Applications	IHC-P
Reactivity	Human

Package

500 µl, 100 µl

PRODUCT

This antibody was validated on 76 different Normal Tissues by IHC-P.

Summary

[Go to Normal Tissue Gallery](#)

[Go to Cancer Tissue Gallery](#)

Applications

Application Note

*Optimal dilutions/concentrations should be determined by the researcher.

Suggested dilution	Recommended dilution
IHC-P	1:100-1:200

Note : Heat-induced antigen retrieval for 5 minutes in an autoclave at 121°C in Antigen retrieval buffer, pH 7.8 (GTX443621).

Not tested in other applications.

Product Note

Highly recommended for IHC-P in human tissues.

Properties

Form	Liquid
Buffer	PBS, 2% BSA, 20% Glycerol
Preservative	0.025% ProClin 300
Storage	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
Concentration	0.25 mg/ml
Immunogen	Recombinant protein encompassing a sequence within the center region of human NFIX. The exact sequence is proprietary.
Purification	Affinity purified by Protein A.
Conjugation	Unconjugated



For full product information, images and publications, please visit our [website](#).


Date 2026 / 10 / 08 Page 1 of 2

For laboratory research use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.

Note

Purchasers shall not, and agree not to enable third parties to, analyze, copy, reverse engineer or otherwise attempt to determine the structure or sequence of the product.


DATA IMAGES

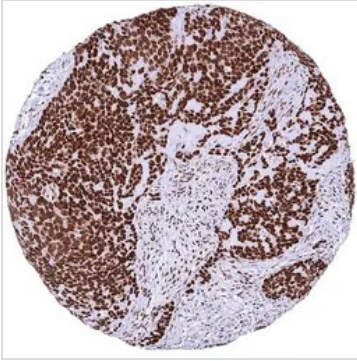

GTX640561 IHC-P Image

IHC-P analysis of human anal canal squamous cell carcinoma (SCC) tissue using GTX640561 NFIX antibody [HMOV329] HistoMAX™.

Squamous cell carcinoma with strong nuclear NFIX staining of tumor cells.

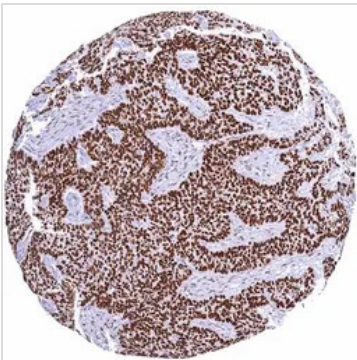

GTX640561 IHC-P Image

IHC-P analysis of human basal cell carcinoma (BCC) from skin tissue using GTX640561 NFIX antibody [HMOV329] HistoMAX™.

Basal cell carcinoma with strong NFIX immunostaining of tumor cells.



GTX640561 IHC-P Image

IHC-P analysis of human fallopian tube mucosa tissue using GTX640561 NFIX antibody [HMOV329] HistoMAX™.

Strong nuclear NFIX staining of all cells.



For full product information, images and publications, please visit our [website](https://www.genetex.com).