

Chk2 (phospho Thr387) antibody

Cat. No. GTX79109

Host	Rabbit
Clonality	Polyclonal
Isotype	IgG
Applications	WB, IHC-P
Reactivity	Human

Package 100 μg

Applications

Application Note

*Optimal dilutions/concentrations should be determined by the researcher.

Suggested dilution	Recommended dilution
WB	1:500 - 1:1000
IHC-P	1:50 - 1:100

Not tested in other applications.

Calculated MW 61 kDa. (Note)

Properties	
Form	Liquid
Buffer	PBS, 150mM NaCl, 0.5% BSA, 50% glycerol (Please contact us for BSA-free format)
Preservative	0.02% Sodium azide
Storage	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
Concentration	Batch dependent (Please refer to the vial label for the specific concentration.)
Immunogen	The antiserum was produced against synthesized phosphopeptide derived from human CHK2 around the phosphorylation site of threonine 387 (C-G-Tp-P-T).
Purification	Purified by sequential chromatography on phospho- and non-phospho-peptide affinity columns. From serum
Conjugation	Unconjugated
Note	For laboratory research use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.
	Purchasers shall not, and agree not to enable third parties to, analyze, copy, reverse engineer or otherwise attempt to determine the structure or sequence of the product.




For full product information, images and publications, please visit our website.

Date 2025 / 12 / 30 Page 1 of 2

€ 886-3-6208988 📻 886-3-6208989 🐷 infoasia@genetex.com




DATA IMAGES

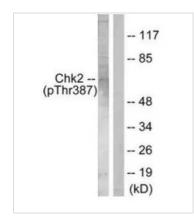

GTX79109 WB Image

WB analysis of Jurkat cell lysate using GTX79109 Chk2 (phospho Thr387) antibody.

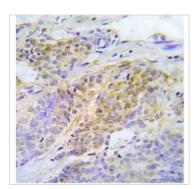

GTX79109 IHC-P Image

IHC-P analysis of human lung adenocarcinoma tissue using GTX79109 Chk2 (phospho Thr387) antibody.



For full product information, images and publications, please visit our <u>website</u>.

Date 2025 / 12 / 30 Page 2 of 2